# SO<sup>®</sup> Wash




oral care system

- EN User guide\*
- IT Manuale d'uso\*
- FR Manuel d'utilisation\*
- DE Bedienungsanleitung\*
- ES Manual de uso\*
- PT Manual de utilização\*
- PL Instrukcja obsługi\*
- RU Руководство пользователя\*
- UA Інструкція для користувача\*
- RO Manual de utilizare\*


\*for "Hydrojet", "Hydropulser" and "Family" version

[www.sowash.it](http://www.sowash.it)

## General Pattern

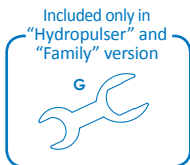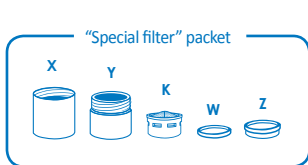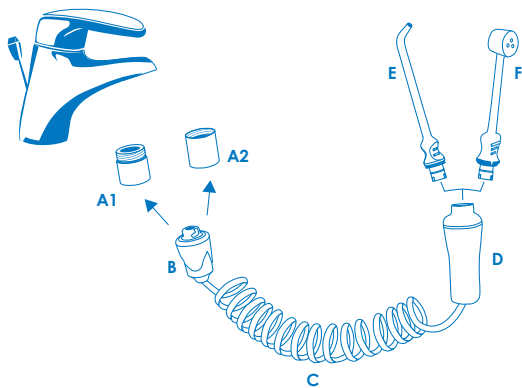
## Filter assembly to the tap

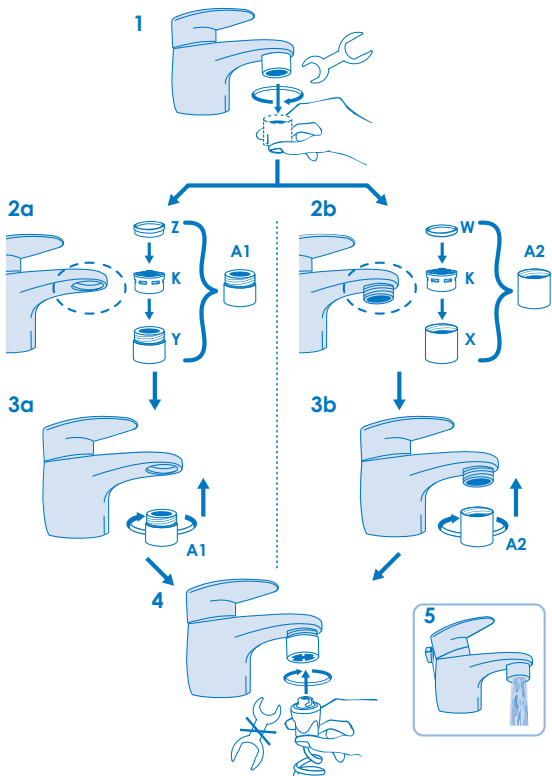Use

6

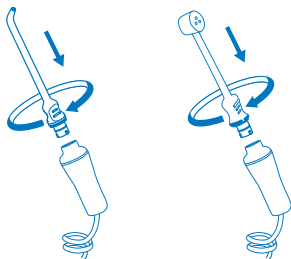
7

8

9

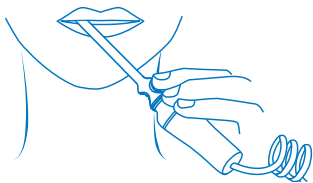
---

10


---

11

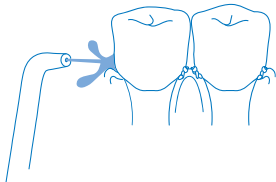5


12

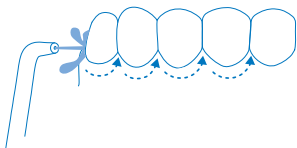
13a

13b

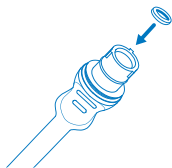14

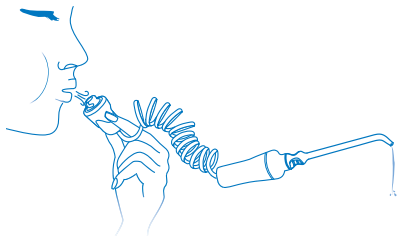**THE PACKAGE CONTAINS:**

- **“Special Filter” Bag containing:**
  - Special filter holder with internal thread (X)
  - Special filter holder with external thread (Y)
  - Aerator Filter (K)
  - Gasket (W - with rectangular section) for Special filter holder with internal thread
  - Gasket (Z - with “step” section) for Special filter holder with external thread
- **Connector (B) with 80cm spiral hose (C) and handle (D)**
- **Hydro-jet Head (E)\***
- **Hydro-pulser Head (F)\*\*/\*\***
- **Key for unscrewing filters holder (G)\*\*\***
- **SoWash carrying bag\*\*\***
- **Multilingual Manual**

*\*1pc in the “Hydro-jet” and “Hydro-pulser” package; 4pcs in “Family” package.*

*\*\*1pc in the “Hydro-pulser” package; 4pcs in “Family” package.*

*\*\*\* included only in the “Hydro-pulser” and “Family” package.*

**FILTER HOLDER ASSEMBLY**

To connect the SoWash aerator filter to the tap (bath sink, shower\*\* or bath tub):

- 1 - Unscrew the aerator filter holder of your tap and make sure that the gasket is removed (fig.1). If necessary, use the key (G - included only in the “Hydro-pulser” and “Family” package).
- 2a - In case of taps with internal thread (fig.2a), insert the aerator filter (K) inside the special filter holder with external thread (Y) and place the gasket with “step” section (Z) on the aerator filter (fig.2a). Screw then to the tap the filter holder with external thread (A1) just assembled (fig.3a).
- 2b - In case of taps with external thread (fig.2b), insert the gasket with

rectangular section (W) and the aerator filter (K) inside the special filter holder with internal thread (X) (fig.2b). Screw then to the tap the filter holder with internal thread (A2) just assembled (fig.3b).

- 3 - Screw the new SoWash filter holder with dry hands, using the connector (B) and turning anti-clockwise (fig.4).

\*The installation of the filter holder to the tap should be done only the first time you connect SoWash: **the exclusive SoWash filter holder for tap does not need to be removed and replaced after each use** because it allows normal use of the tap (fig.5), providing ventilation and water filtration. (The special filter holder SoWash connects to all the taps with unscrewed filter holder according UNI-EN 246 standards as almost all the taps available on the market).

\*\* Through the “multifunction shower joint” present in the range of accessories and sold separately.

To view the phases of assembly, watch videos on [www.sowash.it](http://www.sowash.it) in the “video” and “how to use it” sections.

**HOW TO USE SOWASH**

1. Mount your favourite accessory on the handle (insert and rotate – fig. 6).
2. Connect SoWash to tap filter holder (A1 / A2) with the quick connect (B), by inserting and turning counterclockwise (fig.7). For a better grip, hold the quick connector with dry hands.
3. Leaning towards the sink, insert the head into your mouth, making sure that the hole for the exit of the water is directed to the teeth.
4. Open the valve gradually until it reaches the desired temperature and pressure.
5. Rinse your mouth and interdental spaces (see “Suggestions for use”).

SoWash is equipped with a safety valve

for automatic control of the water flow (B - patented system). **The water pressure in excess is automatically released into the basin (fig.8) to ensure an effective and pleasant outgoing jet.** You may hear a sound, like a whistle, when the tap is fully open. To eliminate this sound just decrease the water pressure from the tap control.

#### **SUGGESTIONS FOR USE**

- Before using it, connect SoWash to the tap and let hot water flow down to the basin for 10-15 seconds to rinse the inside of the hose.
- While using SoWash, keep your mouth slightly ajar to prevent splashing and let the water escape, being bent towards the sink (fig.9 - fig. 10).
- Direct the jet perpendicular to the gumline and clean all areas around the teeth and between a tooth and the other (fig.11 - fig.12), both on the front and inside the dental arches.
- After using SoWash, blow inside the quick connector (B) to release the remaining water in the hose (fig.14) and dry SoWash, especially in the case of sporadic use.

#### **WARNINGS**

- We do not recommend the use of SoWash by children under the age of 12 years, except in the presence of an adult.
- Do not use with water systems with pressure over 6 BAR.
- Do not direct the jet of water under the tongue, to the eyes, nose, ears or other sensitive areas.
- Do not use in presence of open sores on the tongue or in the mouth.
- Do not use water too hot or too cold, it is recommended to keep the temperature lukewarm.
- Do not bite on the heads during use, or tear them with your teeth.

- If there is a leakage of water from the back of the handle (fig.13a), check that:
  - 1) the head is well placed on the handle;
  - 2) the little black gasket is inside mounted head (fig.13b).
- At the end of use, close the water flow before releasing the quick connector.
- Use SoWash only as directed in these instructions.

#### **SPARE PARTS AND ACCESSORIES**

You can purchase the entire SoWash product line on [www.sowash.it](http://www.sowash.it).

[Hydro-jet](#) with single water jet.

[Hydro-pulser](#) with triple pulsing jet.

[Universal Adapter "Addy"](#) allows SoWash use without unscrewing the tap filter holder. It is ideal while travelling, in hotels or in case of non-standard tap threads.

[Multi-function joint for shower tube](#) to connect SoWash to the shower tube

[SoProx Interdental brush 3 in 1](#) can be used:

1. as a normal interdental brush;
2. applied on the practical handle;
3. on the SoWash system as a Hydro-interdental brush.

[SoProx Interdental Brush Spare](#) available in different sizes and shapes for every requirement.

[SoWash water-jet toothbrush](#) with medium/hard bristles.

[SoWash Metal Filter holder](#) available with standard thread, internal or external.

[Adapters](#) for taps with non-standard threads.

#### **WARRANTY**

CSP srl guarantees the main body of SoWash (quick connector B + hose C + handle D) for 3 years from the date of purchase. Are explicitly excluded from the warranty adapters and filter holder for the tap, heads and all accessories. To take advantage of the warranty, the customer must present to CSP

srl the original sales receipt as proof of date of purchase and, if required, the product itself. The warranty is not valid if the product has been misused or has been altered or damaged after purchase. CSP srl may also at its discretion replace only parts of the product found to be defective. The installation is the responsibility of the consumer and is not covered by this warranty.

For any clarification or information please visit [www.sowash.it](http://www.sowash.it) or contact support by email at: [customer@cpsrsl.net](mailto:customer@cpsrsl.net)

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CSP S.r.l. - Italy
CSP S.r.l. uses only solar power